

# GUIDE TO USING YOUR BIOTRANSDUCER®

# CRYSTAL WAVE

---



**BIOTRANSDUCER®**

**TECHNOLOGY**


**FOR**

**TENNANT**

**BIOMODULATOR®**

BioTransducer® technology is a new approach to health care for the 21st Century. The BioTransducer® is powered by your Tennant BioModulator®. It provides therapy that addresses the underlying cause of pain and inflammation without drugs or invasive medical procedures. The BioTransducer® is not a laser.

The BioTransducer® can be used through clothing allowing the BioModulator® therapies to be delivered anywhere on the body without the need of direct skin contact.

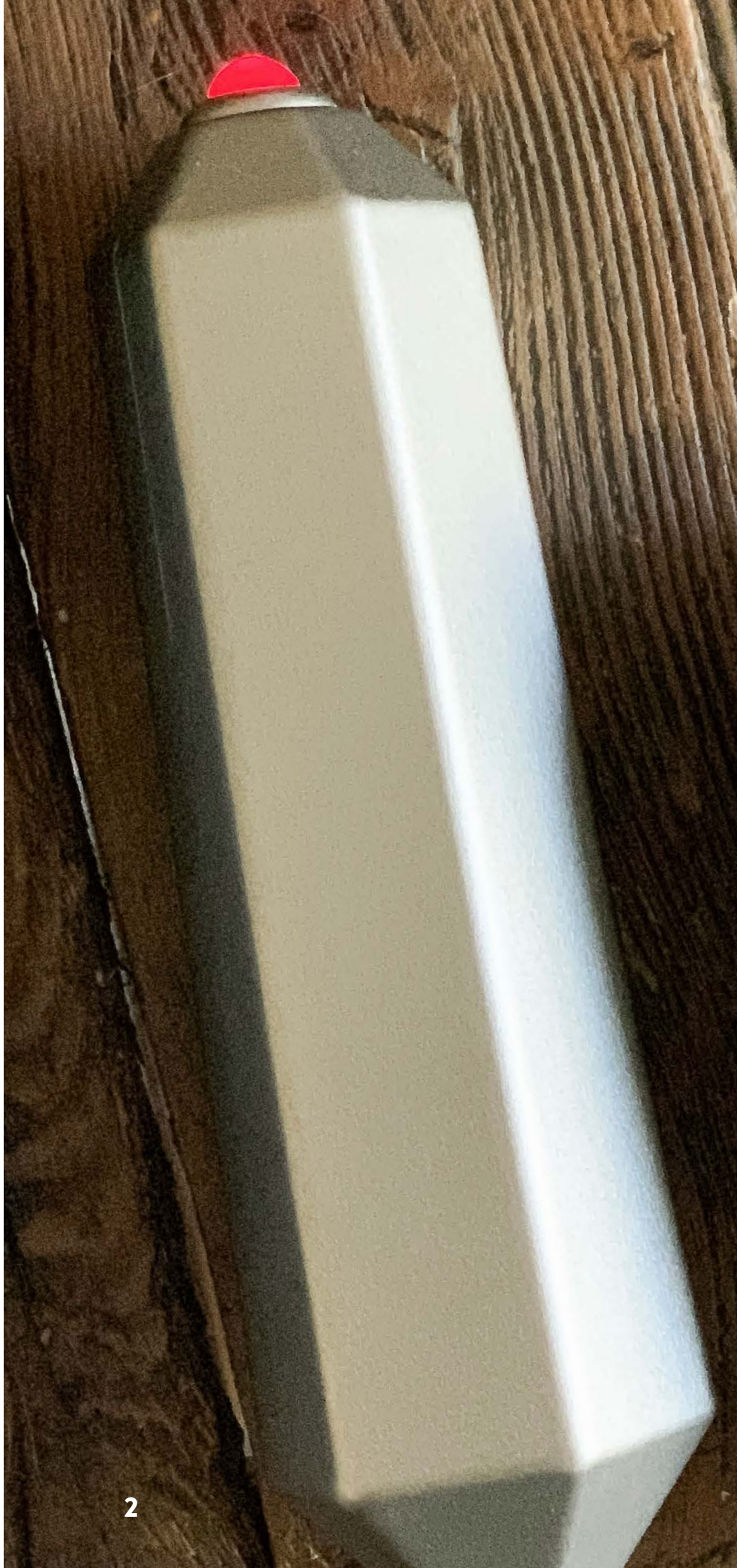

## **BIOTRANSDUCER® THERAPY**

BioTransducer® Technology provides no-touch, painless therapy. The Tennant BioModulator® transmits frequencies and voltage that the BioTransducer® delivers on a carrier wave into the tissue, through skin and bone. Using the BioTransducer® to deliver the frequencies and voltage speeds up therapy time. A therapy that would have been 30 minutes with the BioModulator®.

## **SPECIAL NOTES**

Some users experience a subtle tingle or sensation of warmth from the BioTransducer®. Most people don't feel any sensations when using the BioTransducer®.

You may feel initial apprehension about providing therapy to yourself or others. When you carefully follow the instructions and directions, there's no need to be concerned. The devices and therapy steps were created with the knowledge that application will be general and does not have to be precisely aligned with any given point of area. Even if you are providing the therapy to a different area than intended, it is not harmful and, in most cases, will still benefit you.

# INSTRUCTIONS FOR USE

1. Connect the BioTransducer® to the Tennant BioModulator® with the connection wire.  
*Use care when plugging the wire into the port on the BioModulator®. You can cause damage to your BioModulator® if it is knocked or bumped when plugged in. It is very important to remove the plug when not in use.*
2. Turn on your BioModulator® (Pro or Plus).
3. Set the mode:
  - **For BioModulator® Plus:**
    - **TEN-8** for pain and inflammation related to muscle, bone, ligament and tendons.
    - **INFINITY** for all other issues related to organs and general health.
    - **INFINITY** for eye treatment.
  - **For BioModulator® Pro:**
    - **TEN-8** for pain and inflammation related to muscle, bone, ligament and tendons.
    - **INFINITY** for all other issues related to organs and general health.
    - **VA-STEP** for eye treatment.
4. Adjust the power by pressing the + and - buttons. Since there is no skin contact you may increase the power as high as you want when using it on the body.  
*When using it on the eyes set the power as is indicated in the eye treatment section (pg 6).*
5. Point the BioTransducer® at the area of concern, holding it between 1/2 an inch to 3 inches away from the body.
6. As a general guideline, therapy times are flexible just as they are when using the BioModulator® alone.  
For example :
  - a) 5 – 20 minutes for minor issues and for general use.
  - b) As long as needed for more serious conditions.*Tip: Therapy can continue for several hours when you mount the BioTransducer® to a mic stand and position it on the area of concern.*

*When treating an injury it is important to start treatment immediately.  
Immediate treatments have shown a significant difference in recovery time for patients.*





# EYE THERAPY

# TREATING EYES WITH THE BIOTRANSDUCER® PRO

## SETTINGS:



**FOR PRO**  
Set Mode to  
**VA-Step.**  
Set Power to  
**250.**

OR



**FOR PLUS**  
Set Mode to  
**Infinity.**  
Set Power to  
**solid orange dot.**



## POSITION:

Hold the BioTransducer® Crystal Wave against the temples (as shown).



## DURATION:

Treat **both sides** for **15 minutes each.**



Senergy Medical Group, 9901 Valley Ranch Parkway East, Ste 1009, Irving TX 75063  
(972) 580-0545 • [www.senergy.us](http://www.senergy.us) • [live.well@senergy.us](mailto:live.well@senergy.us)